Female External Catheter

Designed to gently conform to the body.

Indications for Use

The Surity™ Female External Catheter is intended for non-invasive management of urinary incontinence in adults with female anatomy.


Contraindications: Users with urinary retention.

SYMBOL	SYMBOL TITLE	EXPLANATORY TEXT
	Manufacturer	Indicates the medical device manufacturer
	Lot number	Indicates the manufacturer's batch code so that the batch or lot can be identified
	Reference number	Indicates the manufacturer's reference number so that the medical device can be identified
	Medical device	Indicates that this product is a medical device
	Single use	Indicates a medical device that is intended for one single use only
	Do not use if damaged or opened	Indicates that the product should not be used if the package has been damaged or opened
	Use by	Indicates the date after which the medical device is not to be used
	Consult Instructions for Use	Indicates the need for the user to consult the instructions for use


Warnings and Precautions

- The Surity™ Female External Catheter is designed for single-use only. Do not wash or reuse the external catheter. Reuse may result in problems such as skin irritation or infection, malfunction, or product damage.
- Never push or rub the external catheter against the skin during placement or removal to avoid potential skin injury.
- Do not insert the external catheter into the vagina, anus, or other body openings.
- Do not use the external catheter if there is skin irritation or breakdown in the area around the labia or urethra.
- Do not use the external catheter if there are open skin lesions, severe ulceration, or severe rashes in the genital area.
- Do not use during moderate to heavy menstruation if internal menstrual products (e.g., tampons or menstrual cups) cannot be used.
- Do not use if the user has frequent episodes of fecal incontinence.
- Do not use if the user has had recent surgery on the external female anatomy.
- If an allergic reaction occurs, stop use immediately and contact your doctor.

Female Anatomy

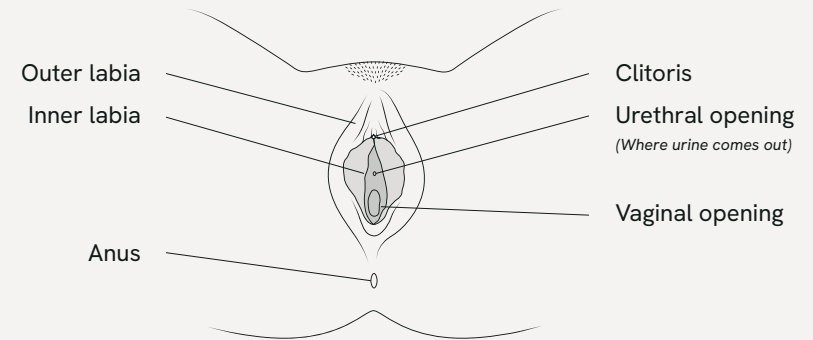Female External Catheter

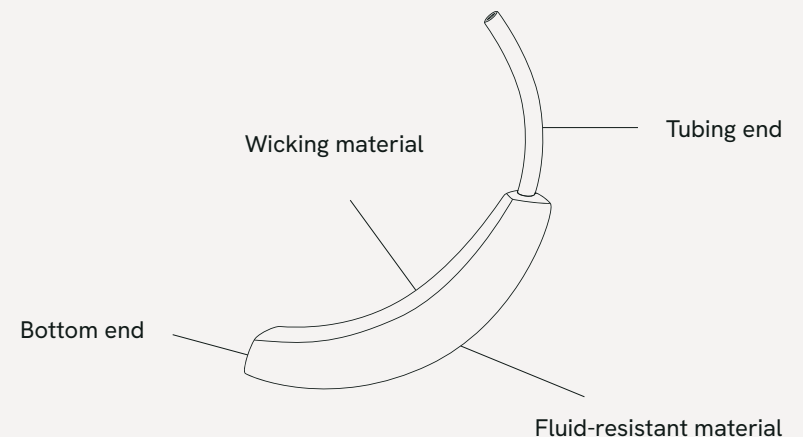Instructions for Use

How to Prepare

- Wash your hands before and after placing or removing the Surity™ Female External Catheter.
- Wear gloves if you're placing or removing this product for someone else.
- Clean the genital area and check the skin condition for rashes and skin breakdown prior to placement.
- Do not use if the user has visible signs of skin irritation or breakdown in the genital area or at the application site.

Urine Collection System Setup

1. Assemble your at-home urine collection system according to its manufacturer instructions.
2. Securely connect the Surity™ Female External Catheter tubing to the collection system tubing, ensuring a snug and leak-free connection. (Fig. A)
3. Turn on the urine collection system by pressing the power switch.

External Catheter Placement

4. Make sure the user is lying on their back or in a reclined position with legs separated.
5. Gently spread the inner labia to expose the urethral opening. Clean the genital area thoroughly from front to back. (Fig. B)
6. Hold the external catheter with the tubing oriented toward the navel (belly button).
7. **Inner labia placement (Recommended):**
Gently tuck the external catheter between the inner labia, placing it directly against the urethra. Make sure the inner labia remain spread and are not tucked underneath the external catheter. (Fig. C)

Outer labia placement (Alternative):

Gently place the catheter on the outside of the outer labia, positioning it against the urethral opening without tucking it into the inner labia—similar to placing a pad.

8. Tuck the bottom end of the external catheter between separated buttocks. (Fig. D)
9. If applicable, gently position underwear or protective briefs over the external catheter.
10. Allow the legs to return to a natural, closed position while ensuring the external catheter stays in place.







External Catheter Removal

11. Separate the legs and gently spread the labia.
12. With the suction still on, gently pull the external catheter away from the body—being careful not to drag it against the skin. (Fig. E)
13. Disconnect the external catheter from the urine collection system tubing. (Fig. F)
14. Dispose of the used external catheter in accordance with local, state, and federal waste disposal guidelines.
15. Turn off the urine collection system by pressing the power switch.

Maintenance and Replacement

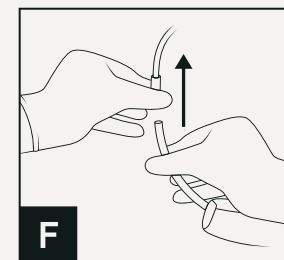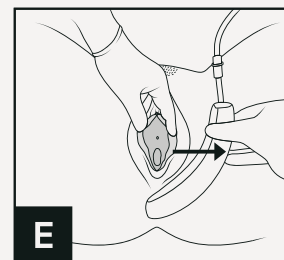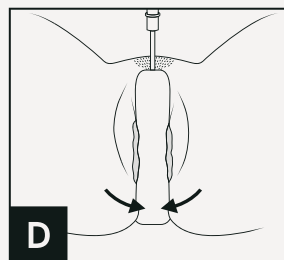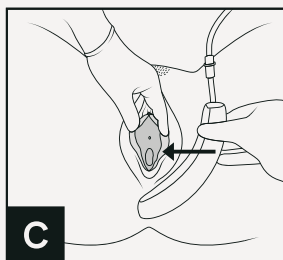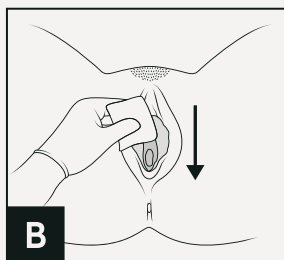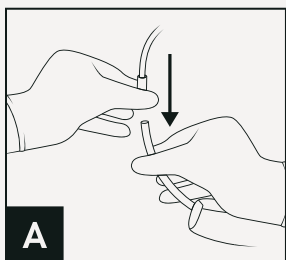- Replace the external catheter at least every 8-12 hours, or immediately if contaminated with feces or blood.
- Always check the skin condition and clean the genital area before placing a new external catheter.
- Store the products in room temperature, 60-77°F.
- If any serious incident arises from the use of this product, contact the manufacturer and national competent authority.

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